




Lymphosep

Lymphosep is designed for a simple, rapid isolation of lymphocytes from whole blood, once diluted and treated with anti-coagulant or defibrinating agent. Lymphosep also allows the rapid isolation of PBMC (peripheral blood mononuclear cells).



- Standardized and high quality medium
- Simple & rapid isolation of lymphocytes from whole blood
- Optimal density of **1.077 g/ml**
- Produced and tested according to the highest quality standards (ISO 13485 and ISO 9001)

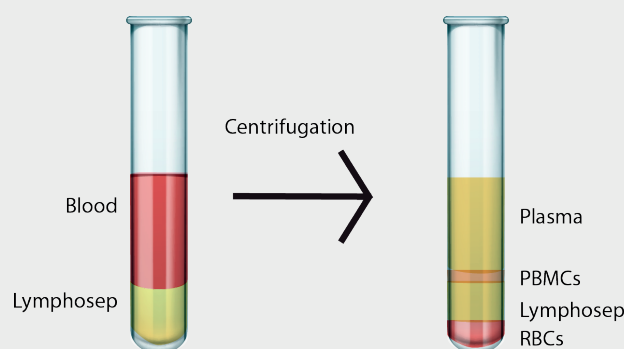

RECOMMANDATION

Use diluted and treated blood (with anticoagulant or defibrinated) for optimum yield.

For best results, use blood drawn less than 2 hours before manipulation.

Do not use with blood drawn 24 hours or more before the manipulation.

Do not use as an API or direct inoculation into human or animals.

Ordering Information

Cat N°	Product Name
L0560-100	Lymphosep, 100 ml
L0560-500	Lymphosep, 500 ml

Please visit www.biowest.net for more information

www.biowest.net

2, rue du Vieux Bourg, 49340 NUAILLÉ

Tel.: (+33) 241 464 242

Email: biowest@biowest.net